



Beyond Scanning : Specialized for All-on-X

FREEDOM **Air**



Meet the IOS That Redefines Possibility

FREEDOM **Air**


Introducing FREEDOM Air, the perfect blend of lightness and precision. It delivers exceptional scan quality, fast scanning speed, and a user-friendly experience designed with hygiene in mind. This all-in-one, complete scanner takes your workflow to the next level in both efficiency and comfort.

Born from countless innovations and tireless development, FREEDOM Air sets a new standard for digital dentistry.

Begin the transformation—with FREEDOM Air.




FREEDOM **Air**

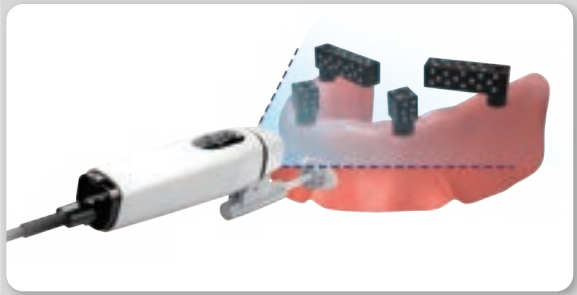

Optimized for All-on-X
Precisely captures implant position data




Smart voice memo via AI recognition
Clinical voice notes tailored to your workflow

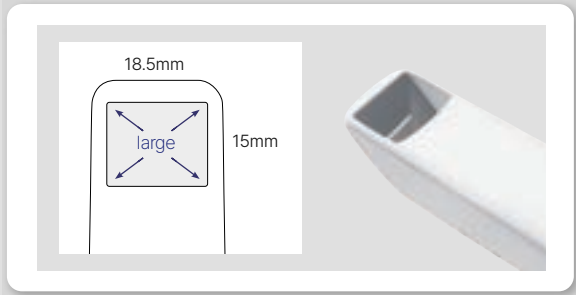

Wide scan field (18.5 × 15mm)
Perfect results in a single scan




Ultra-lightweight at 145g
Comfortable handling with effortless scanning

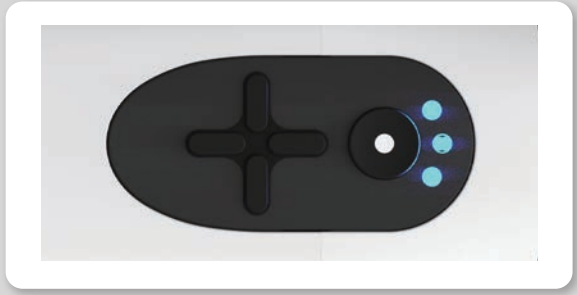

Intuitive D-PAD control
Simple and hygienic operation



LED indicator display
Instantly check scanner status and connectivity

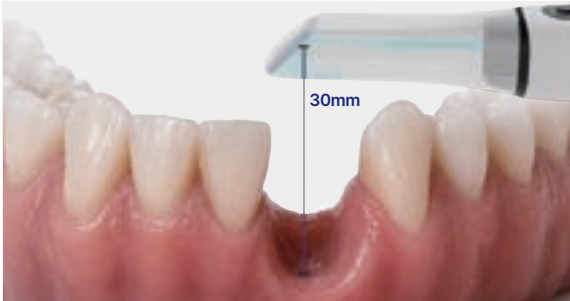

30mm deep scan depth



Realistic oral cavity, ultra-precision real color



Smart scan filter



2-in-1 Scanner

One Device Two Solutions

Extraoral Scan

Photogrammetry

Multi-Implant Position Scan



Intraoral Scan

Stitching Scan

Teeth, Soft tissue Scan



Precision Meets Innovation

A Game Changer : Beyond Scan

PIM Solution

Precise Implant
Measurement Solution

Developed using DOF's proprietary technology, PIM (Precise Implant Measurement) expands the scanning area with a dedicated PIM guide equipped with a special illumination extension lens. This enables simultaneous capture of multiple scan bodies outside the oral cavity.

With wide-area scanning in a single scan, FREEDOM Air captures multiple scan bodies simultaneously without stitching, instantly recognizing their positions for precise implant placement. It is DOF's ultimate All-on-X scanning solution.


Conventional Intraoral Scanners



!

- × Stitching-based Scanning
- × Capturing repeated scans of narrow areas
- × Recognizing the outer shape of the scanbody
- × Higher risk of data distortion

FREEDOM Air PIM Solution

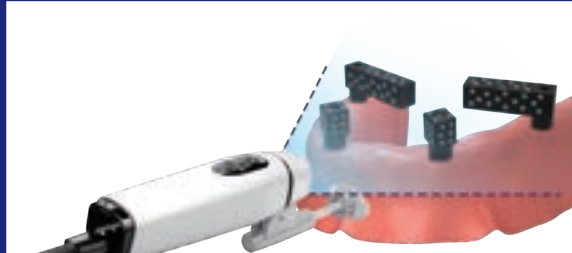

✓

- ✓ Photogrammetry-based Scanning
- ✓ Capturing wide area in a single scan
- ✓ Recognizing the unique coded scanbody
- ✓ Delivering high precision results

Scanning Method

Scan Area


Scanbody Recognition

Accuracy

PIM Only Mode

This mode acquires only precise implant position data through dedicated PIM scanning and aligns it with external intraoral scan data for use.



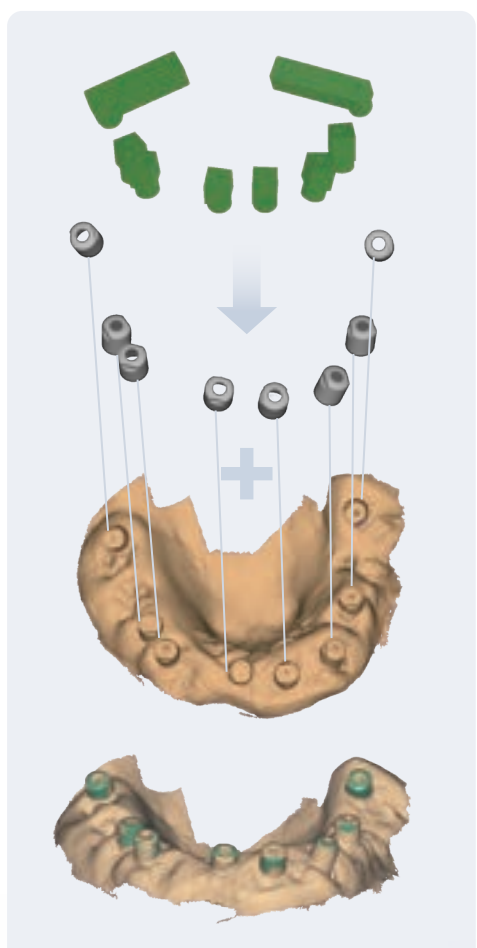

- 1** FREEDOM AIR PIM scan data
→ Convert and extract scanbody geometry data
- 2** Import and align external intraoral scan data with PIM data
- 3** Data alignment completed



Procedure

Ultra-Precise 3D Data Through Photogrammetry

By simultaneously capturing images from multiple angles, FREEDOM Air generates ultra-precise 3D data without distortion. It accurately records the position of scan bodies, ensuring outstanding precision and reliability—even in the most complex full-arch restorations.



PIM Solution

Engineered for Perfect All-on-X



PIM Scanbody

MUA
Multi-Unit Abutment

Screw-integrated scanbody connected to the MUA to capture the implant's 3D position during PIM scanning.

PIM KIT Components




Components		QTY	Image
▪ Scanbody	PG SQ Scanbody (Short, 5*5*9 mm)	6	
	PG L Scanbody (Long, 20*5*9 mm)	2	
	6061 Aluminum (a high-strength aluminum alloy)		
▪ PIM Cap Scanbody	8		
▪ Hand Driver (Hex 1.2)	2	 Short x 1 Long x 1	
▪ Torque Ratchet (35N)	1		

* PIM Scanbody sterilization required before use. (Autoclave at 134°C for 15 min, followed by 30 min drying)

Compatible Multi-Unit Abutment

Type	Compatible MUA
NB	Nobel BioCare
	Straumann (SRA)
	B&B Dental
	Sweden & Martina
	Jdental Care
	DESS
DT	Bio Infinity
	Dentium
	Cowellmedi
DO	C_Tech
NC	DIO
FX	Nucleoss
	FOXDEN

* Type list is continuously updated. * Scan the QR code for the full MUA list.
* For compatibility checks or additional requests, contact DOF customer support.

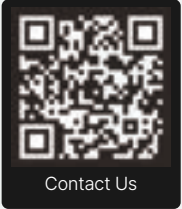


PIM Accessory

Extended PIM Scanbody



An extended scanbody designed for MUAs located below the gingival margin or for cases with swelling after immediate loading. It ensures stable data acquisition even in challenging oral conditions where normal scanbodies are difficult to scan.

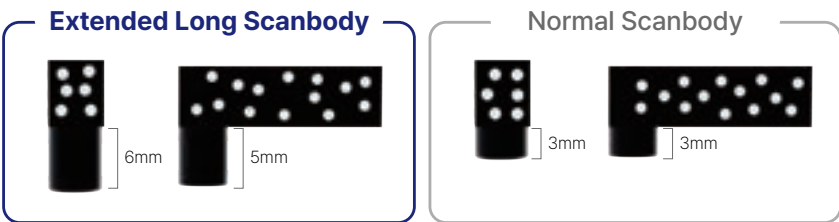

PIM Cap Scanbody



The PIM Cap Scanbody acts as a reference point to align PIM data with soft tissue data. Library-based matching enables fast, highly accurate PIM alignment.



PIM Bite Cap




The PIM Bite Cap mounts to the top of the PIM Scanbody to make bite registration easy. Even for patients without teeth or prostheses, it allows you to capture accurate occlusal data based on the VD (Vertical Dimension).

Size: 3,4,5,6,7mm



PIM Cap Scanbody Scan Workflow

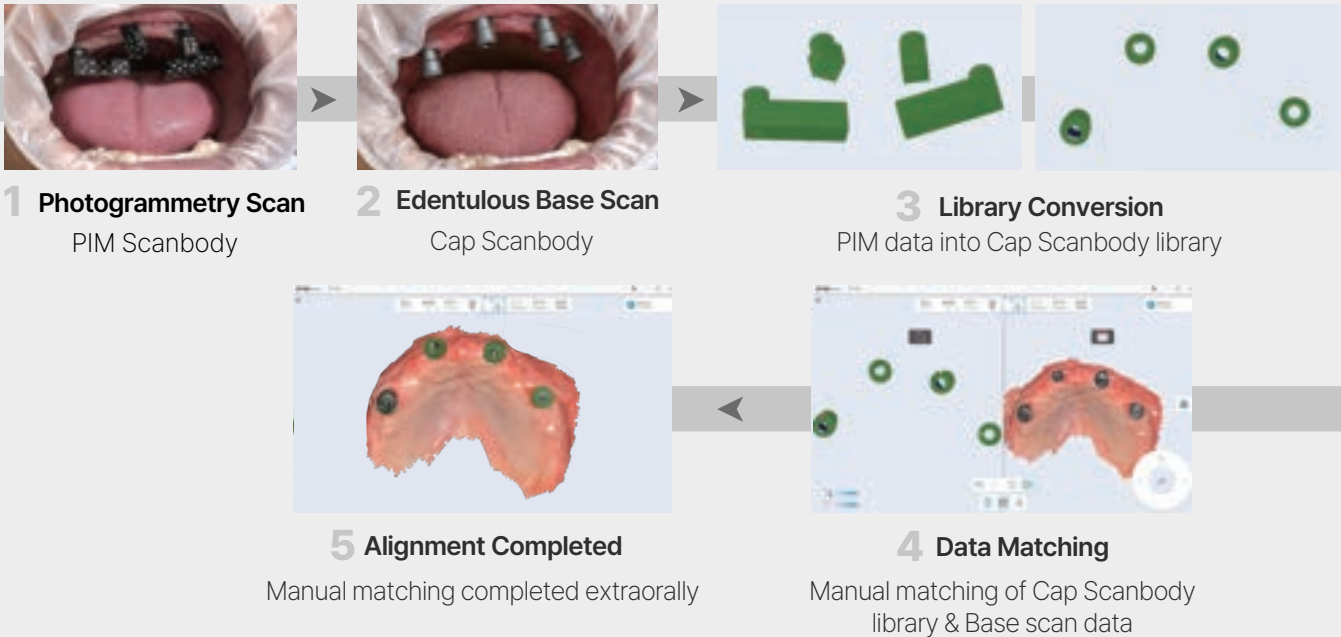



Components

Product Package		PIM KIT Package	
1	FREEDOM Air	1	PIM KIT
2	Small Tip	2	PIM Guide Tool
3	Middle Tip	3	All-on-X Practice Model
4	Large Tip	4	PIM Guide Tool Cover
5	Table Holder		
6	Wall Holder		
7	Cleaning Cotton		
8	Thumb Drive		
9	C to C Cable (1.8m)		
10	A to B Cable (1m)		
11	Power adapter (1.2m)		



Cart Size
W 550 x D 450 x H 840mm
* The Cart and PC are not included and must be purchased separately.

Technical Specification

Product Name		FREEDOM Air
Scanning Technology	Camera Resolution	1920 x 1200 pixels
	Scan Depth	30mm
	Accuracy	10µm
	Anti Fogging	Yes
	Light Source	LED (White, Blue)
	Export File Format	PLY, STL
Scanner Handpiece	Dimensions (Body + Tip)	219 x 35 x 36.6(mm) (W x H x D) without cable
	Weight (Body + Tip)	145g
Connection	Connectivity	USB 3.2 Gen 1
Tip	Tip Size	Middle Tip : 22*17(mm) Large Tip : 25*20.5(mm) Small Tip : 16.3*13.5(mm)
	Field of View	Middle Tip : 16.5*14(mm) Large Tip : 18.5*15(mm) Small Tip : 10*9(mm)
	Autoclavable	Up to 150 times

Recommended System Requirements

FREEDOM Air Recommended System Requirements	
CPU	- Intel 14th Gen i9 or higher - Intel Ultra 7 & Ultra 9
GPU	-NVIDIA GeForce RTX 4060 or higher, RTX 50XX series (VRAM 8GB or higher) -Driver version 576.80 or newer is required (Based on Studio driver)
Memory	Minimum 32GB
Hard Drive	- Minimum 1TB SSD (SSD with DRAM recommended) - At least 20GB of free space
OS	Windows 10/11 (64-bit)

About DOF Inc.

DOF is a CAD/CAM total solution company that develops the world's leading 3D dental scanners and milling machines with proprietary technology. DOF Grows Together with Customers

Founded in 2012, DOF disrupted the industry with its moving-camera scanners and has grown rapidly ever since. With the world's only moving-camera model scanner FREEDOM X5, full-arch optimized intraoral scanner FREEDOM Air, all-in-one milling machine CRAFT 2, compact dry milling machine CRAFT DRY, high-end metal milling

machine CRAFT PRO, and all-in-one titanium-capable CRAFT PLUS, DOF has achieved the true meaning of a model-less total solution. DOF Grows Together with Customers At DOF, every product is designed from the customer's perspective — focusing on required functions and eliminating inconveniences.

Even after release, continuous feedback from users is applied to improve products. DOF's vision and future is to make dental professionals' work more convenient and enjoyable.



doflab.com